

ACCURATE HEMOGLOBIN AND HEMATOCRIT RESULTS WITH ONE SIMPLE TEST

2 results in 1 test

Easy to use

No calibration required

Data management options available



Easy to use


- User friendly features minimises training time
- Step by step instructions on screen
- Backlit touch screen
- User selectable language menu
- Soft-load cuvette holder minimizes risk of contamination

Practical and portable

- Hemoglobin and hematocrit results from one sample available after 25–60 seconds (depending on the concentration)
- 8µl venous, arterial or capillary blood sample
- Microcuvette minimises air bubbles giving consistent results
- Compact in size and weighing just 700g
- No maintenance required. Auto self-test
- Integrated rechargeable battery (100 hours)

Accurate and reliable

- Operating ambient temperature 15°C - 40°C
- Photometric azide methemoglobin method
- Measuring range: 0–25.6 g/dL; 0–15.9 mmol/L
- Imprecision: CV <2 %
- Linearity: 0–20 g/dl: ±0.3 g/dl; >20 g/dl: ±0.7 g/dl
- Calibrated to NCCLS reference method
- Stores up to 4,000 patient results
- Control cuvette included

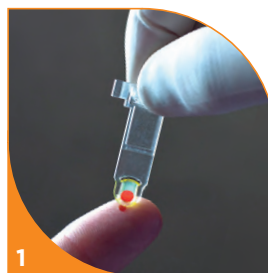

1. Collect 8 µl blood sample.
2. Put microcuvette into analyser
3. Result appears in 25-60 seconds



EKF
DIAGNOSTICS

SCIENCE MADE SIMPLE

SIMPLE SAMPLING FROM ANY ANGLE AND NO AIR BUBBLES

Suitable for
tropical conditions

Convenient package sizes

Optimised shelf life

Low sample volume

Precise sample collection

- Only 8 µl of capillary, venous or arterial blood
- Pre-filled cuvette with dry reagent
- Blood drawn in by capillary action
- Minimised risk of air bubbles give consistent results
- Photometric azide methemoglobin method
- Measuring range: 0–25.6 g/dL; 0–15.9 mmol/L

Individually packaged microcuvettes

- Suitable for tropical conditions
- Carton of 50. Ideal for mobile use and low volume sites
- Each cuvette package contains its own desiccant
- Low sensitivity to humidity
- Operating temperature 15°C–40°C
- Store at room temperature 15–30°C
- Shelf life of 12 months



Distributed by

Manufacturer

EKF diagnostic GmbH
Ebendorfer Chaussee 3
39179 Barleben, Germany
Tel. +49 (0)39203 7850

Sales

EKF Diagnostics Holdings plc
Avon House, 19 Stanwell Road,
Penarth, Cardiff, UK, CF64 2EZ
Tel. +44 (0)29 20 710570
www.ekfdiagnostics.com



Approved for the
European Market

EKF
DIAGNOSTICS

SCIENCE MADE SIMPLE

Disclaimer: EKF-diagnostic GmbH reserves the right
to change the product specifications without prior notice.

Revision 0.0-11/2013